3D smart SHIATSU & TAPPING Massage

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3D smart SHIATSU & TAPPING Massage

User manual

- Thank you for purchasing massage equipment from Axelro AB
- Please keep this manual properly for future reference.
- Before use this appliance read all the instructions carefully.

Product Name: 3D Smart Shiatsu & Tapping Massage

Manufacturer: Axelro AB

Product Model: NMT 0717880238798

Power Output: DC12V 3A

Frequency: 50/60Hz

Timer: 15 min

Weight: 2.2 KG

Size: 42*21,5*22 cm

Thank you for purchasing our 3D Smart SHIATSU & Tapping Massage.

The package includes a massage device, AC adapter, car charger, manual and a protection for the massage balls.

The protection for the massage balls is good to use if you massage directly against the skin so that it does not get so hard, the protection is also washable.

	<p><u>The correct approach for dealing with the product</u></p> <p>This symbol indicates that the product in any EU country is not disposed with other household waste. To avoid any damage to the environment and human health, the product should be recycled. Increase reusability of the remaining available substances and return used appliances. Please follow the return and recycling steps, or contact the local retailers. They can recycle products to be environmental friendly.</p>
--	--

Important Safeguards.

When using electronic products, especially notice and follow the use instructions. Read all instructions before using.

Danger:

To reduce the risk of electrocution, please follow those instructions as below:
Always unplugs or switch off the appliance after using and before cleaning.

Warning:

To reduce the risk of burns, electrocution, fire or injury to persons:

- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
- Don not operates under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.

- Please must be reference for these using instructions to acquaint the purpose of this product when using it. Do not use non-manufacturer offer spare parts for this product.
- Never operate this appliance if has a damaged cord or plug, if it not working properly, if it has been dropped or damaged, or dropped into water. Return it to service centre for examination and repair.
- Do not pull the cord to move the product, use the power cord as a grip.
- Keep cord away from heated surfaces.
- Do not operate the product under the condition of outlet obstruction. Please kep the uotlet open, lint, hair and other objects blocking.
- Do not use this device outdoors.
- Do not use the machine where aerosols (spray bottles) are used or where oxygen is stored.
- Do not allow children or persons with disabilities to use the product without the guidance of a person responsible for safety.
- Take care with the heat function if you have sensitive skin.

Safety regulations.

Cleaning and maintenance of the product.

Before cleaning, make sure the power switch is off and the plug is unplugged. Use a dry cloth to clean the surface of the product.

NOTE!

Do not use abrasive cleaners such as benzene or thinner to clean the product, the color of the product may fade.

Maintenance

Make sure the power switch is off and the power plug is unplugged outlet. Keep the product clean and away from high temperatures.

If the product is not used for a long time, it is recommended to obtain a dust cover to prevent dust from entering and destroying the product.

Keep this manual in a place for future use.

This product is for consumer use only.

Important safety precautions.

Think of.	When moving the product from cold to warm place, please wait 1 hour before using it again to prevent condensation from damaging the machine.
	15 minutes of use is appropriate per occasion.
	Do not use the product with wet hands.
	Do not place the product outdoors.
	Do not use the product with heat on if you are heat sensitive.
	Do not use the product at room temperatures higher than 40 degrees.
	Do not use the power supply with a rated voltage other than that specified for the product, to avoid fire or damage to the machine.
	Do not smoke while using this product.
	Use the AC adapter that comes with the product.
	To avoid damage or accident caused by electricity, moisture or leakage, do not use the product in a humid environment.
	Do not use the product at the same time as other products that have a similar function.
	If you are using the product and feel uncomfortable, cancel and consult a doctor.
	Do not use this product for therapeutic purposes.
	The power cord must not be pinched.
	Do not disassemble the product.
When disconnecting the product, hold the plug. Do not pull the power cord.	
Prohibition	Do not use accessories that may damage the product.
	Do not fall asleep with the product or while using it.
	Do not use this product when you are intoxicated or physically injured.
	Pay attention to any abnormal conditions (such as smoke, odor, etc.).
	Do not use water or liquid when cleaning the product.
	Do not use strong detergents.
	In case of lightning, do not use the product and disconnect the power cord from the electrical outlet.
Power Plug	When servicing, please unplug the power cord from the wall socket and do not use wet hands as it may cause an electric shock or damage.
	After use or when not in use, disconnect the power cord.
	In the event of a power failure, switch off the power switch immediately and unplug the power cord.

Problem solving.

	Phenomenon	Possible cause	Solution
1.	Sounds crackling when turned on.	The sound heard by the friction from the massage balls during normal operation.	No action is needed.
2.	Do not want to start	The plug or cigarette lighter socket is not connected properly.	Check that the plug or cigarette lighter socket is properly connected.
		Incorrect power adapter.	Make sure it is the correct power adapter, the product only works with the original adapter.

If the above does not help and the problems remain, please contact support for your product at axelro.se or nackmassage.se.

Use of the product.

Step 1: Familiarize yourself with the buttons.

Power button:

Press this button once to turn on the device

The medium range shiatsu speed and high range infrared heat is on by default. One minute shiatsu, one minute shiatsu combined with tapping massage will perform and each function will shift after one minute.

Shiatsu lock - start button:

Press this button to lock - start shiatsu massage. This Function is valid only for shiatsu massage functions not for tapping. Direction can not be changed anytime.

Speed control button (default at medium speed):

Press the button 1st time it goes to low speed;

Press the button 2nd time, it goes to high speed;

It will continue to repeat cycles from medium to low to high speed in this sequence.

Heating button:

Heating only.

Heating function can be controlled independently without massage

Press the heating button 1st time, it goes directly to high-intensity Heating.

Press the heating button 2nd time, it goes to low-intensity.

Press the heating button 3rd time, heating is turned off.

Massage & Heating together

Press the Power button and default mode (interchangeably perform shiatsu function or shiatsu & tapping function per every minute, heating ON at high intensity) is activated.

Then press the heating button 1st time to change heat intensity from high to low.

Press the heating button 2nd time to turn off heating.


Tapping lock-start button:


Press this button to lock - start shiatsu & tapping function. The device will perform tapping & shiatsu all the time.

Device will be switched off automatically after 15 minutes when power button or heating button is switched on.


Step 2: Connect the power adapter and the massager correctly. And then locate the massager on the body adequately.

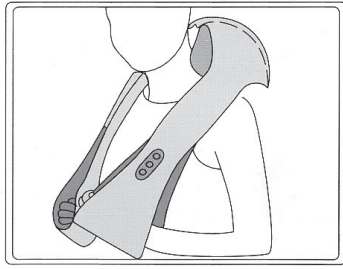
Neck

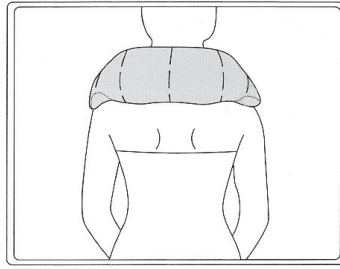
Shoulder

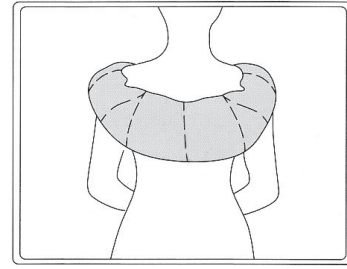
Back

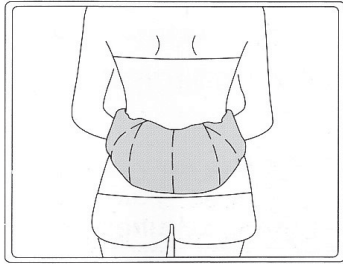Waist

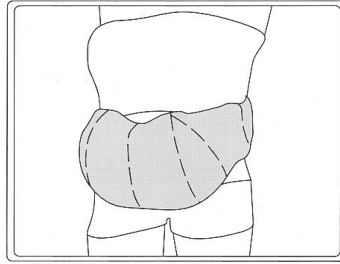Abdomen

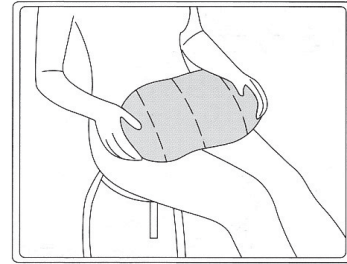Legs

Step 3: Switch on the Power button and enjoy SHIATSU&TAPPING massage. Customize massage functions by selecting Massage direction/Speed/Heating/Tapping buttons respectively.

Note: the unit will be shut-off automatically in 15 minutes.

If you don't want to use the product, please switch off the power button.

Axelro AB
SE559022061101
www.axelro.se
www.nackmassage.se